

Restoration & Prosthetics



Reliable.  
Efficient.  
Precise.  
Alegra.



Turbines

alegra

**TRIPLE SPRAY**

for optimal cooling of the treatment site from all sides.

**SMALL TURBINE HEAD**

for optimal view of the treatment site.

**OPTIMAL LED POSITION**

for optimal view and illumination of the treatment site.

**PUSH-BUTTON CHUCKING SYSTEM**

for simple, rapid bur changes with low actuating force.

**HIGH TORQUE**

for fast removal of tooth structure.

**CERAMIC BALL BEARINGS**

for extremely quiet operation and a long service life.

**PRECISION MANUFACTURE**

for vibration-free and quiet operation.

**WATER FILTER**

for perfect spray function.

**INNOVATIVE GENERATOR SYSTEM**

in the Roto Quick quick coupling for use with Alegria light turbines.




# Long service life and extremely quiet operation

Discover the new Alegra turbines with extremely quiet operation for daily use in your practice. Designed to meet your needs with state-of-the-art technology. Made in Austria: For your best treatment results.



## Long service life

Ceramic ball bearings guarantee minimum wear. Particularly with optimum care and maintenance!

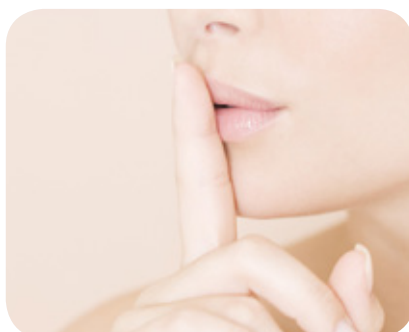

## Quiet

Alegra turbines offer top quality, vibration-free, quiet operation thanks to precision manufacture.



## Replaceable water filter

Prevents dirt particles from entering the turbine.



## Up to the minute

W&H is a dynamic and rapidly expanding company. We guarantee reliable support and outstanding service. Our recipe for success: innovative technologies, high quality material concepts, stringent quality criteria and many years of experience in the development of precision dental instruments.



## Your reliable partner

This combination makes W&H turbines a reliable partner when working on your patients.




W&H Augmented Reality (AR) combines the real world and the virtual world – a unique experience. The future starts now!

## Experience Alegra in a new dimension!

Download the »W&H AR« App to your smartphone or tablet PC free of charge from Google Play (for Android) or the Apple App Store. Open the App and hold your device 50 cm from the picture.

# Maximum efficiency and safety during treatment

The Alegra range of turbines offers high-end technology for top preparation results! Experience maximum efficiency during use combined with maximum safety: for you and your patients.

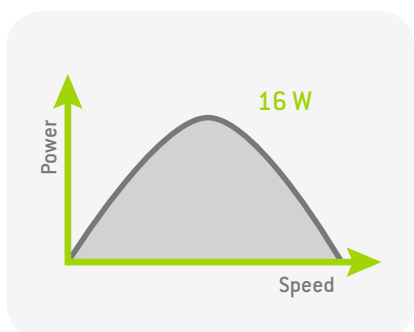


## Fast and precise preparation

Extremely powerful turbines for efficient operation.



## Precision cooling

The triple spray during treatment cools the treatment site from all sides and guarantees safety during treatment.

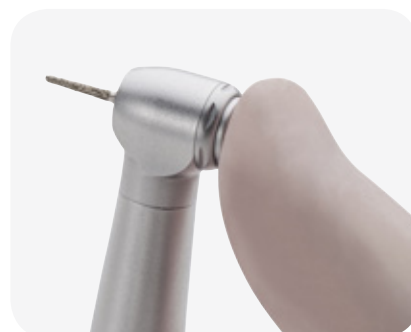

## Simple bur changing

Optimum balance between easy operation of the push-button and maximum retention force of the chucking system. No tools required!



## Tailor-made turbines

From hard-to-reach areas to fast removal of tooth structure: The right head size for every application.



## Patented hygienic head system

Prevents aerosol particles from being sucked into the turbine as the rotor runs down.



# Hygiene & maintenance of your turbines



»The better the care, the longer the pleasure.«  
Regular cleaning and care of your instruments protects them against premature wear and unnecessary service costs.



## For every application

W&H Service Oil F1 is a product of W&H research. It was developed hand in hand with the leading edge W&H drive technology. This oil gives you the security of knowing that you are providing your W&H products with the best possible lubrication.

## Perfect for manual lubrication


W&H Service Oil F1 in the spray can cleans and lubricates your turbines to guarantee troublefree use.



## Maintenance at the press of a button

Assistina 301 plus, Assistina 3x2, Assistina 3x3: The range of W&H cleaning and maintenance units provides the perfect preparation for sterilization from optimum oil service to automatic internal and external cleaning. And all in record time!

## Short lubrication intervals – long service life

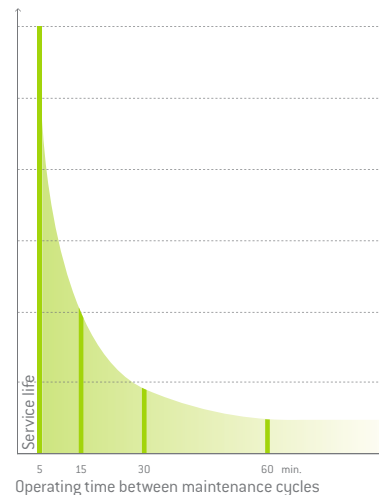

# Optimal view and independent LED

The integration of a generator in the Roto Quick quick coupling gives you a low-cost option of adding an LED turbine to your selection of Alegra turbines. The advantages for you: Flexibility in use and in hygienic preparation.



## Optimal view

The position of the LED in the turbine head guarantees optimal view and illumination of the treatment site.



## LED turbine + quick coupling with generator = independent light

A generator in the Roto Quick quick coupling enables the use of light turbines on the turbine hose without a power supply.



## Excellent ergonomics

Fatigue-free working thanks to easy 360° rotation of the turbine on the coupling, low overall weight and slender design.

## Further advantages at a glance

- > Integrated resuction stop:  
Prevents resuction of aerosol particles into the spray channel
- > Component replacement on site by the user
- > Thermo washer disinfectable
- > Sterilizable up to 135 °C

# Technical data



Scan the QR code with your mobile phone for more information about Alegra turbines

Pictures are for illustrative purposes only. Additional equipment and accessories shown are not included as standard.



Turbines with LED	TE-97 LQ*	TE-98 LQ*	
Couplings with generator	RQ-53 / RQ-54	RQ-53 / RQ-54	
Turbines without light	TE-97 / TE-97 RM / TE-97 BC	TE-98 / TE-98 RM / TE-98 BC	TE-95 RM / TE-95 BC
Light (LED version):	LED and generator	LED and generator	–
Light intensity (LED version):	25,000 lux	25,000 lux	–
Head size:	Ø 10.4 mm	Ø 12.2 mm	Ø 12.2 mm
Power:	14 W	16 W	16 W
Speed:	390,000 min <sup>-1</sup>	330,000 min <sup>-1</sup>	330,000 min <sup>-1</sup>
Rotary instruments:	FG bur, in acc. with ISO 1797-1, Ø 1.6 mm		
Maximum working part diameter:	2 mm	2 mm	2 mm
Maximum length:	21 mm	25 mm	25 mm
Weight:			
> Turbines with LED	46 g	48 g	
> Couplings with generator	26 g / 20 g	26 g / 20 g	
> Turbines without light	46 g / 47 g / 47 g	48 g / 50 g / 50 g	49 g / 49 g
Operating pressure:	2.2 – 2.8 bar (32 – 40 PSI)		
Ball bearings:	Ceramic	Ceramic	Steel
Spray:	Triple	Triple	Single
Volume:	62 dB (A)		

## Your Alegra turbines

Attractively and securely packed.



\* TE-97 LQ and TE-98 LQ turbines with LED are intended for use with W&H RQ-53 and RQ-54 Roto Quick couplings only.

The power and speed of air drives depend on the quality of the turbine hose and may differ from the specified value.

W&H Dentalwerk Bürmoos GmbH  
 Ignaz-Glaser-Straße 53, Postfach 1  
 5111 Bürmoos, Austria  
 t + 43 6274 6236-0  
 f + 43 6274 6236-55  
 office@wh.com  
 wh.com